Digital thermometer instructions

Safety instructions

1. The thermometer should only be used under the supervision of an adult.
2. Do not walk, run or talk during measurements.
3. Store the product safely in the bag when it is not being used.
4. Do not store the product in an area exposed to direct sunlight, dust or humidity. Avoid extreme temperatures.
5. Dropping or subjecting your thermometer to strong shocks should be avoided.
6. Do not attempt to disassemble the unit, except to replace the battery.

Cleaning instructions


Before and after each use, clean the thermometer with a soft cloth and isopropyl alcohol diluted with water, or wash the tip with COLD soap water. DO NOT soak for an extended period of time or sterilize.

What is a normal temperature?

Although the generally accepted “normal” temperature reading is around 36.7°C, temperature readings may vary between 36.1°C and 37.2°C and still may be considered “normal”. Variations in temperature can be attributed to activities such as exercise, smoking, eating and drinking. Your temperature is lower in the morning than in the afternoon. Other variations may be due to the type of measurement. While the oral temperature follows the guidelines previously described, a rectal temperature is generally 0.5°C higher. Conversely, an underarm temperature will be 0.5°C lower.

Use


Features of your digital thermometer

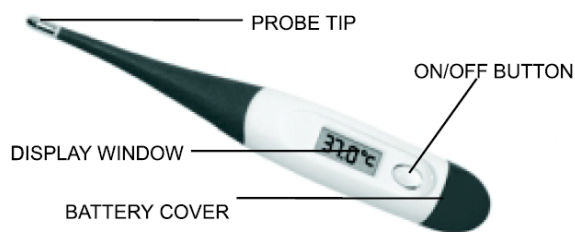1. Press ON/OFF button to activate. The unit will beep and display **188.8°C**. Afterwards, the last measured temperature will be displayed with an “M” mark in the right corner. Then, 37.0°C will appear.
2. The thermometer is ready to start the measurement when **Lo** and a flashing °C are displayed. If the room temperature is higher than 32°C, the room temperature will be displayed instead of **Lo**.
3. When the measurement is complete, you will hear a beeping sound.
4. Press the ON/OFF button to turn the thermometer off. (The unit will automatically power off after approximately 10 minutes.)

Taking an oral temperature

Place the probe tip under the tongue as near as possible to the “heat pocket” (see drawing). Approximate measuring time: 45 seconds.

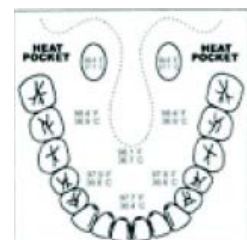Taking an underarm temperature

Wipe the underarm with a dry towel and place the probe tip under the arm so that it is in direct contact with the skin and so that the thermometer is perpendicular to the body. Position the arm across the chest so the probe tip is well covered by the arm. Approximate measuring time: 1 minute

Taking a rectal temperature


Apply a water-soluble lubricant to the probe cover. Gently insert the probe (maximum 2 cm) into the rectum. Approximate measuring time: 40 seconds.

Replacing the battery

The battery in your digital thermometer needs to be replaced when “**T**” appears in the lower right corner of the display. Replace the battery with a 1.55V, LR41 alkaline or SR41 silver oxide or equivalent battery type.

To replace the battery:

1. Remove the battery cover
2. Use a non-metal pointed tool to remove the battery and insert the new battery with the positive (+) side facing up.
3. Carefully slide the battery cover back into place.

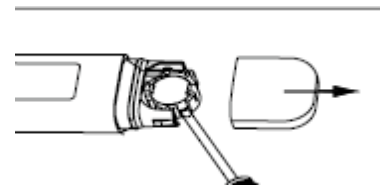ATTENTION: Make sure to keep batteries out of the reach of children and do not store them in hot places!

Specifications:

Measurement range	32.0°C – 42.9°C
Low temperature display < 32.0°C	Lo°C
High temperature display > 42.9°C	Hi°C
Accuracy	35.0°C – 39°C = +0.1°C Other = +0.2°C
Indication of fever (>37.8°C)	30 short beep sounds in 10 seconds
Indication normal temperature (<=37.8°C)	10 long beep sounds in 10 seconds
Auto power off	After about 10 minutes
Battery	LR41 or SR41
Life of battery	2 years
Operating conditions	Temperature: 10°C – 40°C; Humidity: 30% - 85%
Storage conditions	Temperature: 10°C – 60°C; Humidity: 25% - 90%

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